



SYN

Beautiful skin without surgery

ERGY

CACI

SYNERGY+

caci-international.co.uk

THE NEXT GENERATION
IN CACI FACIAL TECHNOLOGY



ADVANCED
NON-SURGICAL
TREATMENTS

Over
30 years
of expertise

ABOUT CACI INTERNATIONAL

CACI International has revolutionised the beauty industry since its beginning in 1992, providing leading edge technology and maintaining an incomparable reputation for excellence and innovation. We've paved the way in anti-ageing, skin rejuvenation and problematic skin solutions for over three decades.

As the UK's number one equipment manufacturer and supplier of non-invasive aesthetic treatment systems, we supply to over 10,000 of the world's leading hotels, spas and aesthetic clinics in 100 countries.

We use the latest technological and scientific advancements to offer a range of treatment systems, which have been specially developed to provide tailored treatments for men and women of all ages and skin types.

Our treatments are often described as the 'red-carpet beauty secret' and have attracted an impressive celebrity following including royalty, film stars, supermodels, beauty editors and make-up artists from around the globe.

THE CACI DIFFERENCE

The multi-functional Synergy is our flagship system and has gained the reputation as the most versatile treatment system on the market. Synergy is our most advanced facial system that has been clinically proven to **reduce wrinkle depth by up to 75% and improve skin elasticity by up to 88%**.

The Synergy features a colour touch screen interface providing a selection menu with 110 treatment programs, including treatments for facial toning, skin exfoliation, skin rejuvenation, problematic skin, sun damaged skin, the reduction of lines and wrinkles, and much more so you can be sure to have a wide variety of treatments for everyone.

The system offers a range of exclusive technologies including our unique **S.P.E.D® Dual Action Technology** that provides simultaneous skin rejuvenation and facial toning.

Synergy is clinically proven* to provide even faster, better and longer lasting results than our earlier CACI facial systems.

CACI is a brand that has become synonymous with delivering effective treatment results without the need for surgery.

When you buy a CACI system you buy more than just a treatment machine, you become part of a winning team and share in the success of the CACI brand name. A CACI system will attract new clients, drive treatment revenues and help your business to grow.

We also help to support your business by providing a wide range of promotional and marketing materials including a marketing guide with tips on how to promote your business, posters, treatment brochures, window stickers as well as a variety of downloadable online resources.

We have a dedicated team of business consultants to advise and help your business to succeed.

At CACI we are committed to being your partner.

*Cutest Systems Ltd 2017 Comparative Microcurrent Study.

SYNERGY+

“The action of two or more things working in unison to produce a combined effect greater than the sum of their separate effects”

CAMBRIDGE ENGLISH DICTIONARY

Whilst CACI systems incorporate multiple technologies, at the heart of all the systems is our core Microcurrent technology.

What sets our Microcurrent technology apart from other Microcurrent machines is our use of ultra-low frequencies in combination with skin resistance technology and a patented digital waveform.


This unique digital waveform optimises results and comfort by gradually increasing in intensity to overcome skin resistance. It then rises to a peak in order to deliver the full power of the Microcurrent energy deep into the muscle. This means that CACI systems not only achieve better results but also longer lasting results.

With the introduction of the Synergy system we have taken Microcurrent technology to the next level.

Synergy's new S.P.E.D® Dual Action Technology represents a significant breakthrough in Microcurrent therapy.

S.P.E.D® enables skincare specialists for the first time to deliver the power of LED light stimulation at the same time as using Microcurrent. This dual action technology has been clinically proven to deliver faster, more effective and longer lasting results than standard Microcurrent.

The introduction of S.P.E.D® means that the best Microcurrent treatment systems on the market are now even better.

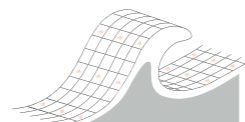

THE CACI
DIFFERENCE

HYDRATONE ROLLERS WITH
S.P.E.D® TECHNOLOGY



MICROCURRENT QUAD PROBES
WITH S.P.E.D® TECHNOLOGY



MICROCURRENT PROBES WITH
S.P.E.D® TECHNOLOGY



ULTRASONIC ACTUATOR



ORBITAL ABRADER WITH
LED LIGHT THERAPY



WRINKLE COMB WITH
LED LIGHT THERAPY



T BAR




110
Treatment
Programs


INCLUDING:
FACIAL TONING
JOWL LIFTING
LINES / WRINKLES
EYE BAGS
SKIN HYDRATION
SKIN PEELING
ACNE / BLEMISHES
STRETCH MARKS
SCARRING
NECK TREATMENTS
CELLULITE
HAND TREATMENTS



FOUR UNIQUE SYNERGISTIC TECHNOLOGIES:

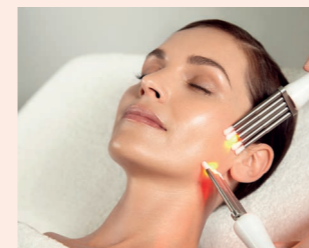

S.P.E.D® Microcurrent
Microcurrent toning
combined with simultaneous
LED light therapy.

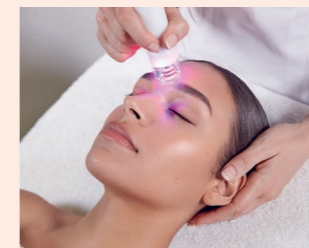

Orbital Skin Resurfacing
Crystal free dermabrasion
combined with simultaneous
LED light therapy.



Wrinkle Reduction
Collagen stimulating waves
combined with simultaneous
LED light therapy.



Ultrasonic Rejuvenation
Piezoelectric dermal
vibrations combined with
simultaneous
ultrasound technology.

WITH STATE OF THE ART DESIGN AND INCORPORATING THE LATEST ADVANCEMENTS IN ELECTRONICS, THE CACI SYNERGY SYSTEM IS AT THE FOREFRONT OF AESTHETIC TREATMENT TECHNOLOGY.

Our continuing investment into research and development enables us to create exciting new and exclusive technologies that deliver dramatic results without the need for surgery or invasive procedures.



Microcurrent Probes with S.P.E.D® Dual Action Technology

The dual-tipped Microcurrent Probes utilise CACI's pioneering S.P.E.D® Dual Action Technology, allowing skincare specialists to harness the power of LED light therapy at the same time as using microcurrent. The dual action of simultaneous LED and microcurrent energy stimulates tissue regeneration, helps in the production of collagen, tightens sagging muscles whilst smoothing the appearance of fine lines and wrinkles. The unique combination of these technologies provides more visible and longer lasting results.



Microcurrent Quad Probes with S.P.E.D® Dual Action Technology


The Microcurrent Quad Probes utilise CACI's pioneering S.P.E.D® Dual Action Technology. They are designed to grip a larger surface area and work the jaw muscles within the Jowl Lift treatment. This treatment specifically targets muscle laxity around the jawline which can cause drooping jowls and loose skin. The treatment has been developed to lift and firm muscles and redefine facial contours to give a firmer, more toned appearance.

Orbital Abrader

This new exfoliation technology uses a handpiece with an abrasive tip that rotates at different speeds on an elliptical axis. The orbital rotation of the abrasive head performs a progressive and controlled skin exfoliation, removing imperfections and impurities. A selection of single use, disposable abrasive treatment tips coated with a layer of silicon carbide are used to either superficially or deeply exfoliate the epidermis, resulting in a clearer, smoother and more youthful looking complexion. Precision control over the depth of abrasion allows you to specifically tailor each client treatment. This unique crystal and vacuum free method of microdermabrasion overcomes the blockage, maintenance and downtime issues commonly associated with crystal based microdermabrasion systems.

A massage tip can also be inserted onto the handpiece to provide a lymphatic drainage treatment.

The Orbital Abrader delivers simultaneous LED light at the same time as the exfoliation to promote healing and skin rejuvenation.

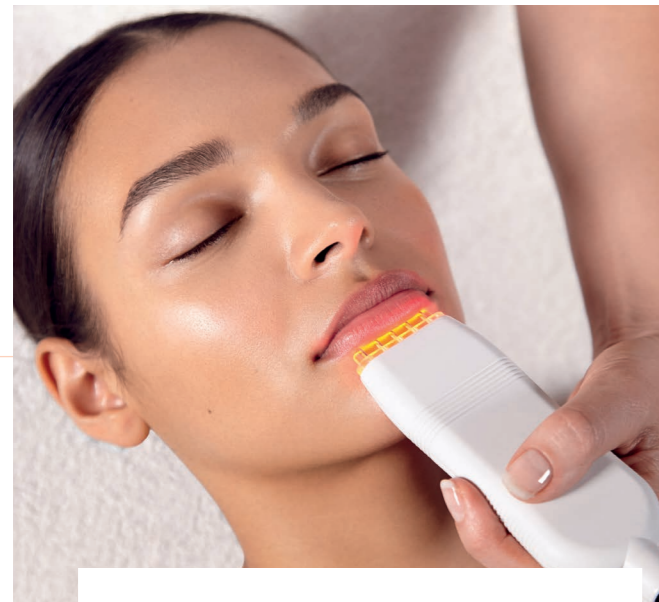

Wrinkle Comb

The Wrinkle Comb is used to effectively and painlessly plump out and soften deep lines and wrinkles, thereby providing a needle free alternative to dermal fillers.


The unique handset transmits high frequency electrical stimulation simultaneously with LED light, which stimulates an increase in elastin and the production of collagen. Red LED light encourages skin rejuvenation and promotes tissue healing. Blue LED Light has an antibacterial and calming effect on inflamed or irritated skin. An alternating pulsed LED mode is particularly effective at treating acne prone or problematic skin.



Ultrasonic Actuator

The Ultrasonic Actuator handpiece contains ceramic/quartz transducers that when excited by an electrical current emit rapid vibrations, 27,000 per second. The vibrations propel impurities to the skin surface and cause dead skin cells to become dislodged. The handpiece is then used to deeply cleanse and gently exfoliate the epidermis. The Ultrasonic Actuator provides a far gentler and less aggressive method of skin exfoliation than traditional microdermabrasion.

The handpiece can also be set to a specific frequency to enhance the absorption of skin care products using sonophoresis (ultrasound waves), to push the products deeper into the skin.

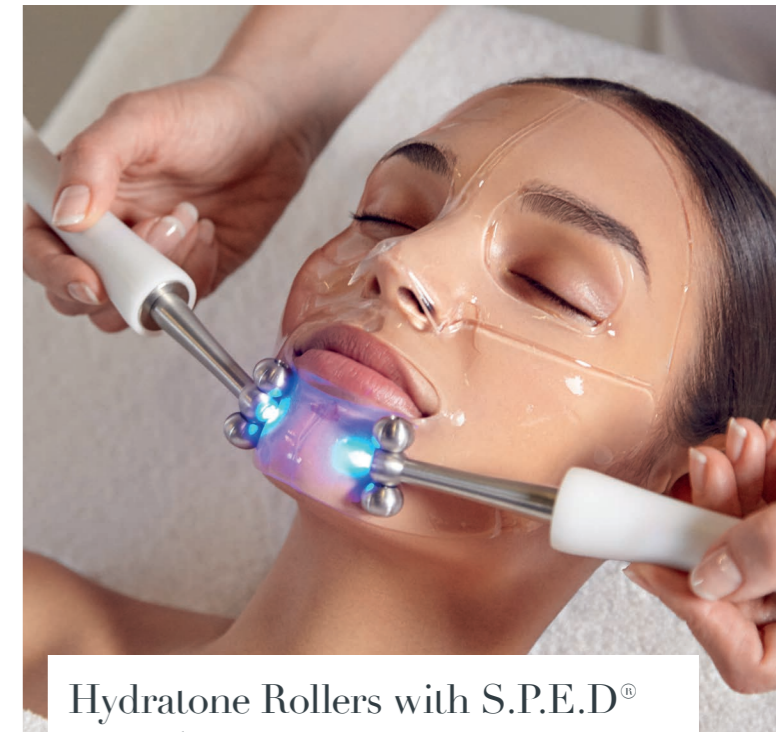

Hydratone Rollers with S.P.E.D® Dual Action Technology

Hydratone S.P.E.D® Dual Action Technology Rollers are used to gently massage the face over the electrically conductive CACI Hydro Mask. The microcurrent's rejuvenating energy transmitted through the rollers provides immediate and intensive skin hydration and helps to restore muscle tone and firmness. Red light is generally used for anti-ageing as it stimulates the production of collagen and rejuvenates the skin. Blue light has a calming action on inflamed or irritated skin.




Enhanced Tissue Repair (ETR) Pads

The Enhanced Tissue Repair face pads are placed on either the face or neck during a non-surgical facial toning treatment, enabling you to work the face & neck at the same time.

A continuous flow of microcurrent is transmitted through the pads, which helps to heal and repair blemishes, scars, stretch marks and pitted skin.

T-Bar

The T-Bar works in conjunction with the Microcurrent Probes to grip a larger muscle surface area, to provide dramatic firming and lifting results on the face.

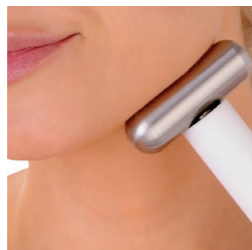

Conductivity Boosters Instant Firming Serum

The Conductivity Boosters are individually filled with a pre-measured amount of Instant Firming Serum specially developed to boost and enhance microcurrent electrical conductivity. The Instant Firming Serum contains active Hexapeptides designed to instantly firm, tone and redefine facial contours. Use these with the Microcurrent Probes to enhance your treatment results.



OPTIONAL ADVANCED ATTACHMENTS



Eye Revive Rollers

The Eye Revive Rollers are designed for treating the delicate eye area. Eye Revive Gel can be placed in the roller hand-piece, allowing the skin to better absorb the product.

The Eye Revive Rollers use microcurrent to gently tighten and tone sagging muscles together with a deeply nourishing CACI Eye Revive Mask to lift hooded eye lids, reduce puffiness and dark circles and soften the appearance of fine lines and wrinkles around the eyes.



Electro Cellulite Massager (ECM)

The ECM combines the therapeutic benefits of massage with the effectiveness of microcurrent therapy. The conductive rollers improve the appearance of cellulite by gently breaking down fatty deposits, kneading tired shoulders, soothing aching limbs and firming sagging muscles. The massaging action of the ECM stimulates circulation and lymphatic drainage, whilst at the same time toning muscles and reducing inches.




Electro Gloves

The Electro Gloves are made of an electrically conductive silver fabric, that when activated stimulate the hands with microcurrent energy.

The CACI Rejuvenating Hand Mask is used in conjunction with the Electro Gloves to nourish and smooth skin texture, improve skin firmness and boost circulation.



WORLD LEADERS IN AESTHETIC TECHNOLOGY

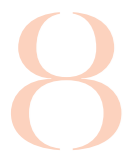




BEAUTIFUL SKIN
WITHOUT SURGERY

RESULTS


Independent Clinical Study

CACI Synergy, introduces pioneering S.P.E.D® Dual Action Technology that delivers simultaneous LED skin rejuvenation and microcurrent facial toning. An independent clinical study conducted by Cutest Systems Ltd, compared CACI Synergy S.P.E.D® enhanced microcurrent with standard microcurrent facial toning. 25 female participants aged between 35-60 took part in the study, 10 treatments were administered with measurements and assessments taken after treatment 1, 5, 10 and 2 weeks following the final treatment.

Skin Elasticity Improvement

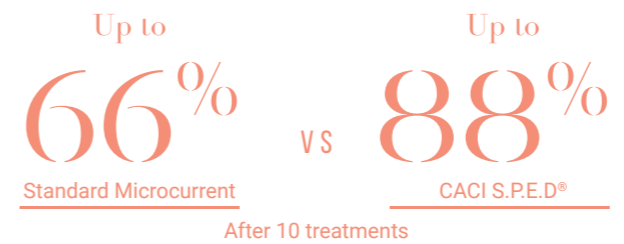

Wrinkle Depth Reduction

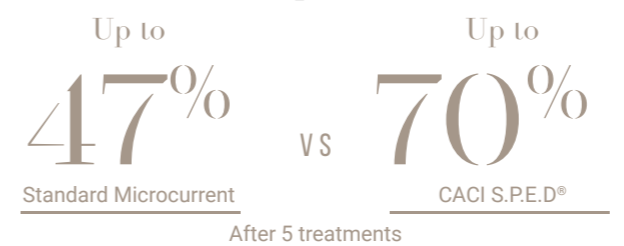

Study Participant Questionnaire



SUMMARY


CACI S.P.E.D® produces faster, better and longer lasting results in improving skin elasticity compared to standard microcurrent.

CACI S.P.E.D® produces faster, better and longer lasting results in the improvement of lines and wrinkles compared to standard microcurrent.


CACI S.P.E.D® results in a faster improvement in skin texture than standard microcurrent.

CACI S.P.E.D® reduces the appearance of blemishes on the skin more effectively than standard microcurrent.


CACI S.P.E.D® provides a greater reduction in red and brown spots than standard microcurrent.

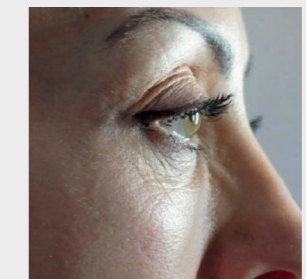

BEFORE

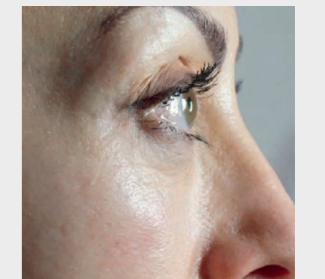

AFTER

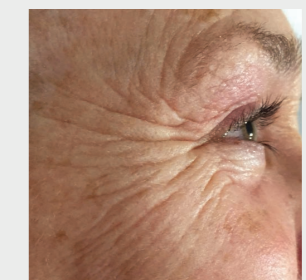

BEFORE



AFTER



BEFORE



AFTER

INCREASE YOUR REVENUE

With Treatments That Appeal to All Clients

Choosing a CACI Synergy gives you the comfort and peace of mind knowing that you are purchasing direct from the UK's largest aesthetic equipment manufacturer and from a brand that you can trust to deliver quality, reliability and performance.

Synergy offers a wide range of treatment options that will appeal to men and women of all ages and skin types.

Offering CACI treatments within your business will enable you to attract new clients, increase client retention and significantly boost your treatment revenues.


With a variety of attractive payment schemes available you don't even need to make a capital outlay and can start generating profit as soon as your CACI Synergy arrives.

Our payment schemes are so readily affordable that you only need to perform just one CACI client treatment course per month to be in profit.

A CACI Synergy really is a must have for any beauty business serious about generating substantial additional income.



SHARE IN THE SUCCESS OF THE CACI BRAND NAME AND INCREASE YOUR REVENUE

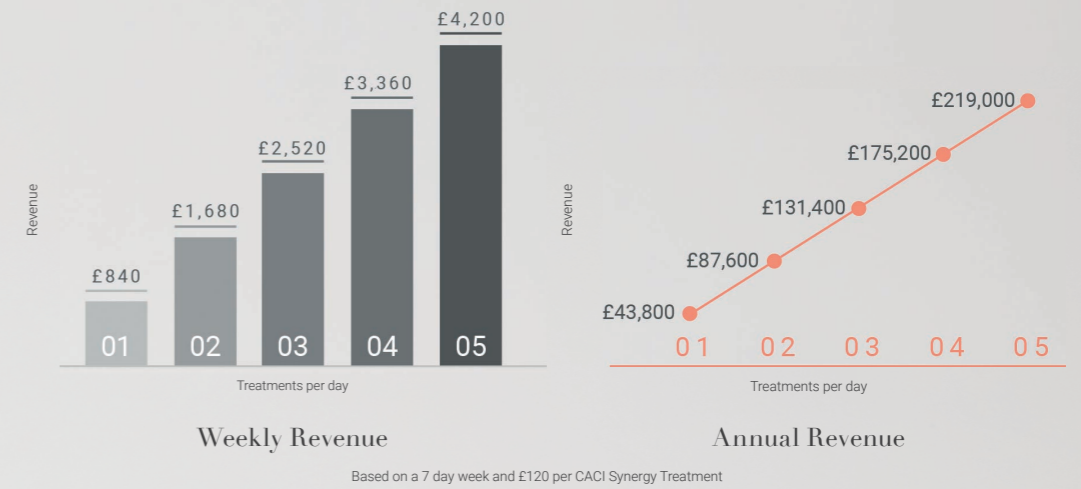

SYNERGY TREATMENTS

TREATMENTS

CACI Signature Non-Surgical Facial Toning (60 Minutes)

This non-invasive facial features S.P.E.D® Dual Action Technology that provides simultaneous skin rejuvenation and facial toning. This treatment uses microcurrent impulses, to lift and tone and LED light therapy to improve skin elasticity and reduce the appearance of fine lines and wrinkles. Firms the face and gives the skin a more youthful, fresh appearance and glow.

CACI Synergy Advanced Non-Surgical Facial Toning (90 Minutes)

Combines CACI's Signature Non-Surgical Facial Toning with advanced skin exfoliation techniques to revitalise the skin, resulting in a brighter and smoother complexion. A combination of advanced technologies are used to plump fine lines and wrinkles, instantly lift and tone the face and firm the neck.

CACI Synergy Purifying (45 Minutes)

This treatment tackles problematic skin and leaves the complexion looking smooth and radiantly healthy. A combination of red and blue LED light therapy, deep cleansing and skin exfoliation techniques are used to effectively brighten and even the skin tone.

CACI Gentleman's Facial (45 Minutes)

Using a combination of advanced technologies, this treatment is designed to specifically address male skincare concerns. Deep cleansing and skin exfoliation techniques tackle razor bumps and in-growing hairs, electrical impulses tighten, firm and redefine the jawline; the CACI Hydro Mask and LED light therapy are then used to hydrate, calm and smooth the skin, leaving the complexion refreshed, revitalised and more youthful looking.

CACI Skin Calm (30 Minutes)

Targets skin prone to breakouts and irritation. This treatment uses LED light therapy that is renowned for its anti-bacterial and anti-inflammatory healing properties, combined with a gel mask to soothe, calm, reduce redness and rehydrate the skin.

EXPRESS TREATMENTS

CACI Hydratone (15 Minutes)

This treatment is designed to deeply nourish, hydrate and firm the skin. Facial rollers emitting simultaneous microcurrent impulses and LED light therapy are gently massaged over a hydrating gel facial mask. The rejuvenating energy transmitted through the rollers provides immediate and intensive hydration and helps to restore skin firmness, soften fine lines and rejuvenate the skin. Ideal for the treatment of dehydrated or sun damaged skin.

CACI Wrinkle Revolution (30 Minutes)

Targets deep lines and wrinkles, such as crow's feet and laughter lines. This treatment is a non-invasive, needle free alternative to collagen injections and dermal fillers that will instantly plump fine lines and wrinkles from the very first treatment.

CACI Lip Booster (15 Minutes)

Enhance your clients' pout. Designed to give fuller, plumper lips and soften top lip lines. This quick fix treatment will have your clients' pout party ready in no time.

CACI Skin Confidence (30 Minutes)

Designed to soften and refine the appearance of stretch marks and scarring, this treatment uses advanced skin exfoliation techniques together with the regenerating effects of microcurrent stimulation and LED light therapy to trigger skin tissue repair and healing.

ADVANCED TREATMENTS

CACI Jowl Lift (15 Minutes)

Target muscle laxity around the jawline using S.P.E.D® Dual Action Technology, which harnesses the power of LED light therapy at the same time as using CACI microcurrent. This treatment will lift and firm the muscles and redefine facial contours to give a firmer, more toned appearance. Ideal for women wanting to improve the appearance of sagging jowls and for men who desire a more chiselled jawline.

CACI Eye Revive (30 Minutes)*

Lift and firm muscles around the eye area, particularly where you may tend to lose definition as you age and develop 'hooded eye lids'. This treatment uses serum filled microcurrent rollers together with the soothing and calming Eye Revive Mask to reduce puffiness and dark circles, and soften the appearance of fine lines and wrinkles around the eyes.

CACI Hand Rejuvenation* (30 Minutes)

The anti-ageing and nourishing hand treatment uses conductive microcurrent gloves with CACI Rejuvenating Hand Mask to reduce fine lines and intensively moisturise, leaving the hands hydrated with a softer more youthful appearance.

CACI Electro Cellulite Massage* (30 Minutes)

Gives an instant lift, shapes and tones the buttock and thigh area. This treatment combines the therapeutic benefits of a massage with the effectiveness of microcurrent therapy. The conductive rollers improve the appearance of cellulite by gently breaking down fatty deposits. It stimulates circulation and lymphatic drainage, which helps to flush toxins away giving a smoother, tighter and dimple free appearance.

Please note: Advanced & Express Treatments can be performed as either stand alone or as add ons to any of our other treatments.

*These treatments can only be performed using the optional advanced attachments.



Technical specification



Dimensions(WDH) : 45cm x 42cm x 15cm

Weight : 9Kg

Touch Screen : 10.4in

Supply Voltage : 60–240 volts

Output Voltage : 30 Volts max (loaded)

Consumption : 110 vA

Safety Class 1 type BF : in accordance with IEC 601 and IEC 601-2-10

Channels : 4 independent

MICROCURRENT WITH S.P.E.D®

Polarity : positive, negative & alternating

Range : 20uA to 640 uA

Frequency : 0.1Hz–800Hz

Waveshapes : 8

Mode : Automatic Sequential

LED BLUE : 475nm – 60mW

LED RED : 640nm – 60mW

WRINKLE COMB

Output : 0.1mA–8mA

Output Frequency : 6000Hz

LEDs : 5x Blue 430nm–40mW

LEDs : 5x Red 640nm–30mW

ORBITAL ABRASION

Rotation Speed : 0–4000 rpm

Motor Voltage : 6–15 Volts

Amperage : 0.1–0.5 Amps

LEDs : 6x Blue 430nm–40mW

LEDs : 6x Red 640nm–30mW

ULTRASONIC HANDSET

Mechanical Oscillation : 30micro +/- 20%

Frequency : 27,000Hz +/-10%

